

Daon pen H5

SPECIFICATION

Brand: Daon pen

Model: H5

Power model: Wireless

Battery capacity: 700 mAh

EMS intensity level: 4 levels

EMS frequency: 0.25 Hz

Speed Level: 6 Speeds

Loadable liquid capacity: 3 ml

Adjustable liquid output: 3 levels

Adjustable needle length: 0 ~ 1.0 mm

Available adapter: US/ EU/ UK/ AU/ KR

Available needle cartridge: H12, H36, H24, Nano-HS, Nano-HR



3 ml




Needle cap regulates liquid output



Uniform Fluid

ADVANTAGES

- Integrates triple effects: EMS electroporation, automatic serum apply and microneedling.
- Portable wireless design, convenient and easy to use.
- Adjustable liquid output 1-3 levels, can adopt essence with different concentration.
- Low battery warning and charging status display L/F, it's more convenient to use.
- Compact one-hand device for simple and efficient operation.
- Ultra-long working time, big capacity battery to ensure the pen working efficiently.

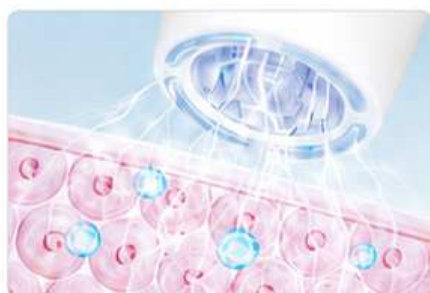

EMS annular micro current
+
carrier essence



One-hand operation



Adjustable needle length 0~1.0 mm

Dr.pen H6

SPECIFICATION

Power model: Wireless

Speed: 6 levels

LED light: Blue, Red

Handle Material: ABS

Battery capacity: 1000 mAh

Adapter: US/ EU/ UK/ AU/ KR


Adjustable needle length: 0~2.0 mm

Liquid Loading Capacity: 3 ml

Automatic Liquid Dispensing Adjustment: 3 Levels

Available needle cartridge: 18 pins, 36 pins,

42pins, Nano-HS, Nano-HR

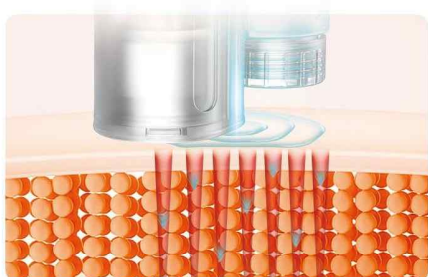

Automatic Penetration



One-Hand Operation



Smart Needle Adjustment



ADVANTAGES

- It is equipped with a liquid-carrying needle and pre-filled with essence before use. During the application process, it automatically penetrates the skin to achieve efficient skin care.
- The flow rate of the essence is precisely controlled, with three adjustable flow rates available to precisely control the treatment time
- It is perfectly suitable for all kinds of thicknesses of essence on the market
- One-click electric needle adjustment eliminates the need for manual mechanical needle adjustment with both hands. The operation process is simple and efficient.



Hydra. Needle 30

SPECIFICATION

Brand: Hydra. Needle

Model: HN30

Needle number: 30 pins

Capacity of vial: 10 ml

Adjustable liquid output: 3 levels

Vial material: Glass

Needle material: Titanium

Available needle size: 0.25 mm, 0.5 mm,
1.0 mm, 1.5 mm

ADVANTAGES

- To meet different thickness of the serum on the market, liquid output is adjustable from level 1 to level 3.
- Equipped with a cap to protect the needles, which keeps needle clean and hygienic.
- The needle head can be purchased separately. It can match standard ampoule on the market. Sealing well and easy to carry.
- Very comfortable to hold with soft silicone, suitable for beauty salon and home use.
- Titanium alloy needles, fine and sharp needles are used for minimal invasive skin care with remarkable results.
- Serum will come out automatically by manual pressing, easy and safe operation for home use/salon/clinic.



Adjustable Liquid Output



Compatible with Ampoules



Deep into the skin



Professional care



Facial care



Hair care



Solve skin problems

ATTIJOU Revitalizing Ampoule_NAD is a 4th-generation Bio-Rebuilder that goes beyond simple nourishment to address the root cause of skin aging: "Cellular Energy Depletion."

By combining clinically proven **PN (Polynucleotide)** for tissue regeneration with **Nicotinamide Adenine Dinucleotide (NAD+)** and **Nicotinamide Mononucleotide (NMN)**, a key driver of cellular energy metabolism, this solution restores the compromised skin barrier and stimulates collagen production within the dermis.

Mechanism of Action

- Mitochondrial Activation:** **NAD+** reactivates mitochondrial function, which declines with age, thereby restoring the cell's intrinsic recovery potential.
- ECM Reconstruction:** The synergistic action of **PN** and **HA** improves the **Extracellular Matrix** (ECM) environment to enhance skin density and elasticity.

Key Benefits (Intended Use)

- Skin exhibiting loss of elasticity and fine lines.
- Aging skin suffering from chronic dryness and lack of vitality.
- Skin with a damaged barrier due to external stimuli.
- Patients seeking overall complexion improvement and a healthy glow.



Treatment Guide

While following standard protocols for premium skin boosters (MTS, Mesotherapy), this guide focuses on maximizing the absorption rate of the NAD component.

Method/Area	Detailed Protocol
Injection (Mesotherapy)	Target Layer: Superficial to Mid Dermis Volume: 0.02~0.05mL per point (Create Micro-papule) Total Volume: Approx. 2.5mL per session (Full Face) Needle: 32~34G (To minimize pain and ensure precision)
Combination Therapy	MTS: Apply immediately after creating micro-channels to induce intradermal absorption. Energy Device: Apply after RF/HIFU procedures for soothing and regeneration boosting.
Treatment Interval	Initial Phase: Every 2 weeks, 2~3 sessions. Maintenance: Every 2~3 months.
Treatment Areas	Face, Neck, Hand

※ Number of treatments and the quantity administered depend on the degree of skin aging.

ATTIJOU Brightening Ampoule

DAON
MEDICAL

ATTIJOU Brightening Ampoule is a high-performance brightening solution that inhibits the entire process of melanogenesis.

It features an optimal blend of **Glutathione** (key ingredient in premium skin brightening IV therapies) and **Tranexamic Acid** (the gold standard for melasma treatment). This formula blocks pigmentation while the **Hyaluronic Acid** base provides an immediate hydration glow, transforming dull and uneven skin tones into a clear, transparent complexion.

Mechanism of Action

- **Melanogenesis Inhibition:** **Tranexamic Acid** and **α -Arbutin** inhibit tyrosinase activity, blocking melanin production at the source.
- **Antioxidant Defense:** **Glutathione**, a potent antioxidant, neutralizes free radicals and protects the skin from oxidative stress.
- **Hydration Glow:** **HA** locks moisture within the dermis, delivering a radiant "water-glow" effect immediately post-procedure.

Key Benefits (Intended Use)

- Skin with hyperpigmentation issues such as melasma and dark spots.
- Dull skin with uneven tone.
- Devitalized skin caused by dehydration.



Treatment Guide

We recommend the **"Dual-Depth Technique"**, adjusting injection depth according to the pigmentation pathology.

Treatment Area	Target & Depth	Interval & Volume
Full Face (Tone-up)	Superficial Dermis • For overall tone improvement & hydration	• Interval: 1.0cm • Volume: 0.01~0.03mL per point
Spot Target (Melasma/Pigment)	Mid-dermis • Intensive targeting to inhibit melanocyte activity	• Interval: 0.5cm(Dense application) • Volume: 0.02mL per point
Tools	Nano-needle or Automatic Injector is recommended for uniform delivery	
Treatment Interval	• Initial Phase: Every 2 weeks, 4~5 sessions • Maintenance: Every 3~6 months	

※ Number of treatments and the quantity administered depend on the degree of skin aging.

The **OLLA Refine** is a specialized Bio-Stimulator designed for delicate skin rejuvenation.

Formulated with **Ultra-Low Molecular Weight PLLA microparticles**, it focuses on improving skin texture and restoring fine details for delicate contouring.

It holds a **total mass of 92mg**, including 50mg of pure PLLA, offering a soft, natural collagen remodeling experience suitable for thin skin layers.

Mechanism of Action

- **Neocollagenesis:** ULMW-PLLA microparticles stimulate fibroblasts to generate new collagen, thickening the dermis and improving elasticity.
- **Dermal Saturation:** The 7mL high-volume format allows for broad, even distribution (saturating) to treat fine lines and pore size across the entire face.
- **Rapid Bio-Integration:** ULMW-PLLA microparticles hydrate quickly and integrate smoothly into the tissue, minimizing downtime and nodule formation.

Key Benefits (Intended Use)

- Patients with thin skin or loss of facial elasticity.
- Skin seeking natural improvement in texture, pores, and fine lines.
- Individuals requiring overall facial rejuvenation without excessive volume bulk.



Treatment Guide

We recommend standard Mesotherapy or precise injection techniques to maximize the dispersion of the fine PLLA particles.

Method/Area	Detailed Protocol
Preparation	Hydration: Inject 5~6mL of Sterile Water for Injection (SWFI) or Saline into the vial. Allow to hydrate for at least 30 minutes. Suspension: Shake vigorously immediately before use to ensure a homogeneous suspension.
Injection (Mesotherapy)	Target Layer: Deep Dermis to Subcutaneous junction. Technique: Micro-droplet, Multi-puncture, or Linear Threading. Cannula (25G–27G): Highly recommended for the smooth delivery of the high-content PLLA suspension. The use of a cannula is the primary method suggested to minimize the risk of vascular occlusion and ensure even distribution of the product. Needle (30G): Optionally applicable for precise, localized treatment. When using a 30G needle, the use of a 'Thin-wall' type is recommended to prevent potential clogging of the suspension
Volume	Point Vol: 0.02 ~ 0.05mL per point. Total Vol: 1 Vial (Full Face or targeted areas like neck/hands).
Treatment Interval	Initial Phase: Every 3~4 weeks, 3 sessions. Maintenance: Every 6~12 months.

※ Number of treatments and the quantity administered depend on the degree of skin aging.

※ Professional Use Only: This product must be administered exclusively by trained medical professionals who are thoroughly familiar with PLLA injection techniques and facial anatomy.

The OLLA Forte is a high-performance collagen stimulator formulated for intensive volumizing and deep structural repair.

With a **total mass of 276mg** and a high concentration of 150mg pure PLLA, it delivers three times the collagen-stimulating power of the Refine version, effectively addressing deep wrinkles and significant volume loss.

Mechanism of Action

- **Structural Scaffolding:** The high-density ULMW-PLLA microparticles create an immediate scaffold in the subcutaneous layer, supporting collapsed tissue structure.
- **Intensive Volumizing:** Induces robust collagen production over time, filling deep hollows and correcting facial contours.
- **Long-Lasting Restoration:** The particles degrade slowly and predictably, maintaining volume and elasticity for an extended period.

Key Benefits (Intended Use)

- Deep wrinkles such as nasolabial folds and marionette lines.
- Areas with significant volume loss (sunken cheeks, temples).
- Facial contouring and lifting for sagging skin.



Treatment Guide

Due to the high concentration of particles, the "Dual-Depth Technique" or use of a cannula is recommended to ensure safety and uniformity.

Treatment Area	Target & Depth	Interval & Volume
Preparation	Hydration: Inject 5~6mL of Sterile Water for Injection (SWFI) or Saline into the vial. Allow to hydrate for at least 30 minutes. Suspension: Shake vigorously immediately before use to ensure a homogeneous suspension.	
Deep Wrinkles & Folds	Subcutaneous Layer Structural support and deep filling.	Volume: 0.1~0.2mL per point (Retrograde injection) Tool: 25G~27G Cannula recommended
Global Volumizing (Cheeks, Temples)	Deep Subcutaneous / Periosteum Volume restoration.	Volume: 1 Vial Technique: Fanning or Cross-hatching.
Tools	Cannula (25G-27G) is strongly recommended to prevent vascular occlusion and ensure smooth delivery of the high-content suspension.	
Treatment Interval	Initial Phase: Every 4~6 weeks, 2~3 sessions. Maintenance: Every 12~18 months.	

※ Number of treatments and the quantity administered depend on the degree of skin aging.

※ Professional Use Only: This product must be administered exclusively by trained medical professionals who are thoroughly familiar with PLLA injection techniques and facial anatomy.



DAON
M E D I C A L

DAON MEDICAL CO.,LTD.

M_ +82-10-7565-5737

E brian@daonmedical.co.kr

www.daonmedical.co.kr

Sales Office #1013, LG Eclat 2nd Seohyun-dong, Bundang_gu,
Seognam-si, Gyeonggi-do, Korea 136591

Factory #1303~#1305, Daeryung Technotown 17th, 25 Gasan Digital
1-ro, Geumcheon-gu, Seoul, Korea 08594